


GSM60-V2T4 Datasheet

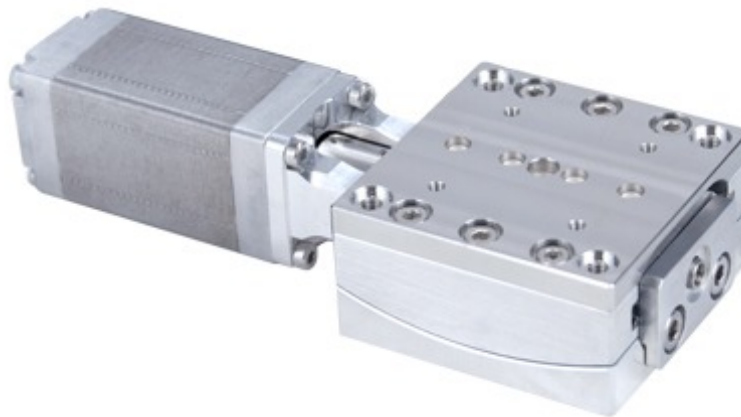

- Vacuum compatible to 10^{-6} Torr
- 40 mm and 60 mm centre of rotation options allow for two axes of motion about a common centre of rotation
- Crossed roller bearings for smooth motion and high stiffness
- Designed for use with an X-MCC Series stepper motor controller or any 2-phase stepper motor controller
- Only 7 feedthrough wires required to control these units

GSM-V2 Series Overview

For more information about the basics of a vacuum system and considerations to keep in mind when gathering requirements for your application, read our technical article, "Motion Device Design Considerations for Vacuum Applications".

Zaber's GSM-V2 Series devices are high-vacuum compatible, computer-controlled, motorized goniometers with a compact size.


Like all of Zaber's products, the GSM-V2 Series is designed to be very easy to set up and operate. If you are considering a multi-axis system, in the XY configuration, these stages are designed to stack with mounting alignment features in the stage top and base.

These stages are wired with a male D-sub 15 connector for plug-and-play use with our X-MCC Series

stepper motor controllers. Simply cut the connector off for flying leads in vacuum and match the wire connections outside of the vacuum for an easy connection. Our handy kits include free software, and all of the accessories that you will need to get the stage running right out of the box.

For more information visit: <https://www.zaber.com/products/vacuum-compatible-stages/GSM-V2>

GSM-V2 Series Part Numbering & Options

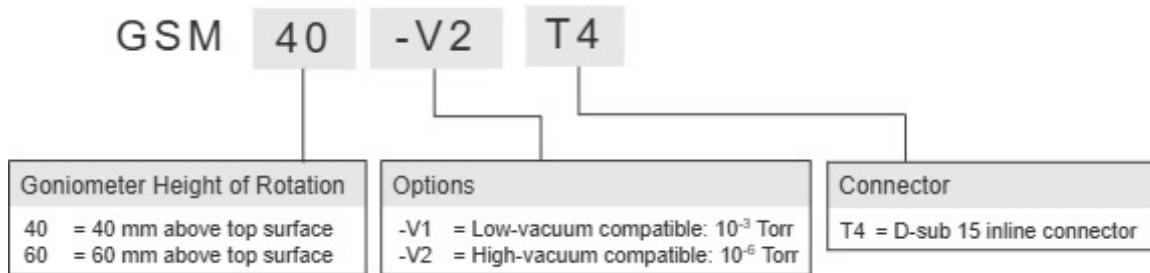

GSM60-V2T4 Drawings

- [GSM-V.png \(Drawing for the GSM-V\)](#)


GSM60-V2T4 Specifications

Microstep Size (Default Resolution)	0.0001131413171° (1.975 µrad)
Built-in Controller	No
Recommended Controller	MCC (48 V) Recommended
AutoDetect	No
Range	+/-7°
Accuracy (unidirectional)	0.06° (1.047000 mrad)
Repeatability	< 0.005° (< 0.087 mrad)
Backlash	< 0.01° (< 0.174 mrad)
Maximum Speed	30°/s (5 rpm)
Minimum Speed	0.000069°/s (1.204 µrad/s)
Speed Resolution	0.000069°/s (1.204 µrad/s)
Encoder Type	No
Maximum Centered Load	30 N (6.7 lb)
Rotation Centre Height	60 mm (2.362")
Rotation Centre Accuracy	< 0.1 mm (< 0.004")
Angular Motion Per Motor Rev	1.44820885888°
Motor Steps Per Rev	200
Motor Type	Stepper (2 phase)
Motor Rated Current	800 mA/phase
Inductance	3.5 mH/phase
Motor Connection	Teflon flying leads with D-sub15 connector
Default Resolution	1/64 of a step
Motor Frame Size	NEMA 08
Guide Type	Crossed roller bearing
Mechanical Drive System	Precision Worm Gear
Limit or Home Sensing	Magnetic home and away sensors
Manual Control	No
Axes of Motion	1
Mounting Interface	M2 and M3 mounting holes
Operating Temperature Range	0-50 °C
CE Compliant	Yes

Microstep Size (Default Resolution)	0.0001131413171° (1.975 μrad)
Vacuum Compatible	High vacuum (10 ⁻⁶ Torr)
Weight	0.320 kg (0.705 lb)

GSM-V2 Series Charts

Recommended Duty Cycle

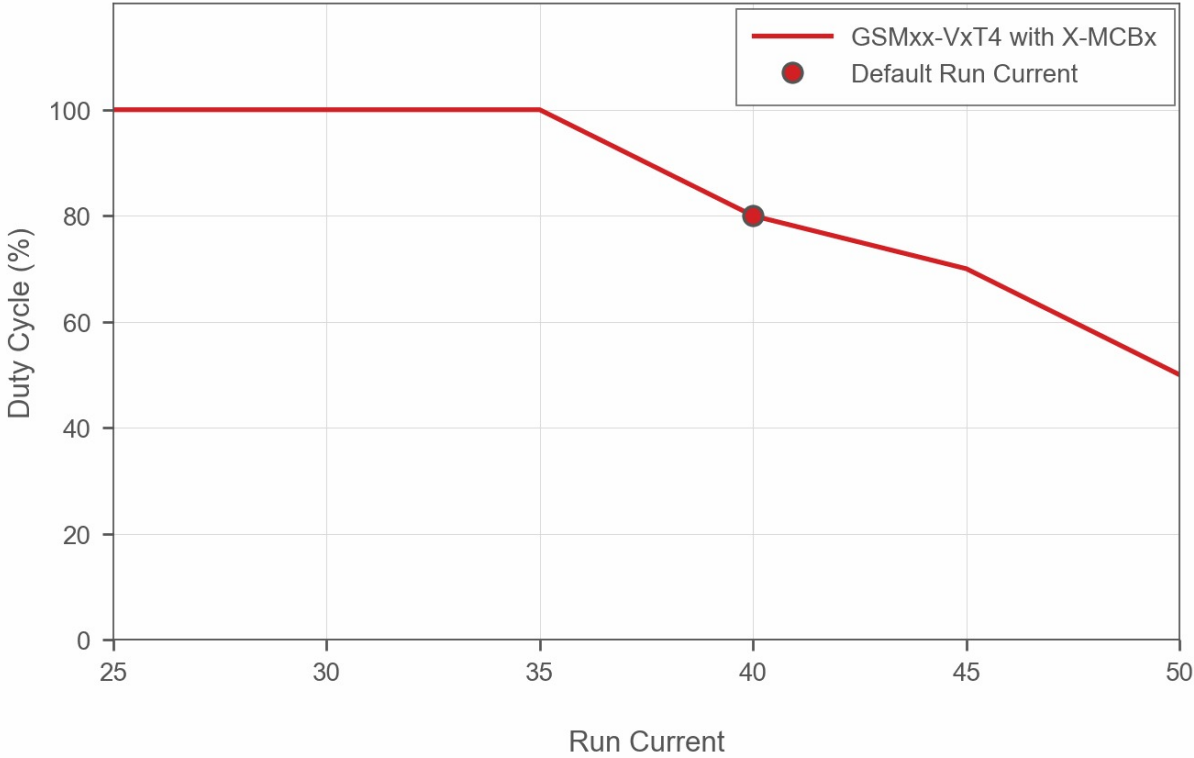

Contact

Email: contact@zaber.com

Phone (toll free Canada/USA): 1-888-276-8033

Phone (direct): 1-604-569-3780

Fax: 1-604-648-8033

Zaber Technologies Inc.

#2 - 605 West Kent Ave. N.

Vancouver, British Columbia

Canada, V6P 6T7

<https://www.zaber.com>